

Entecavir

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IE [4.6 x 250 mm, 5 µm]

Part No: 85325

Mobile phase: Ethanol : Methanol : DEA (90/10/0.1, v/v/v)


Flowrate: 0.7 mL/min

Column Oven Temperature: 25 °C

UV Detection: 260 nm

Diluent: Mobile phase

Sample Concentration: 0.25 mg/mL

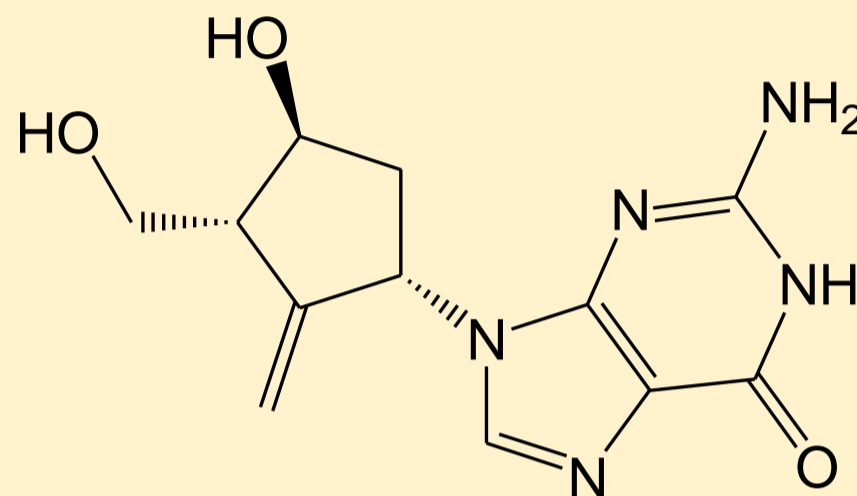

SEPARATION PERFORMANCE

t_{R1} : 5.8	N_1 : 4648	T_1 : 1.5	
t_{R2} : 7.7	N_2 : 3324	T_2 : 1.7	R_s : 4.2

