

Tetrabenazine

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IC-3 [4.6 x 250 mm, 3 µm]

Part No: 83525

Mobile phase: Methanol : ACN : DEA (90/10/0.1, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 280 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

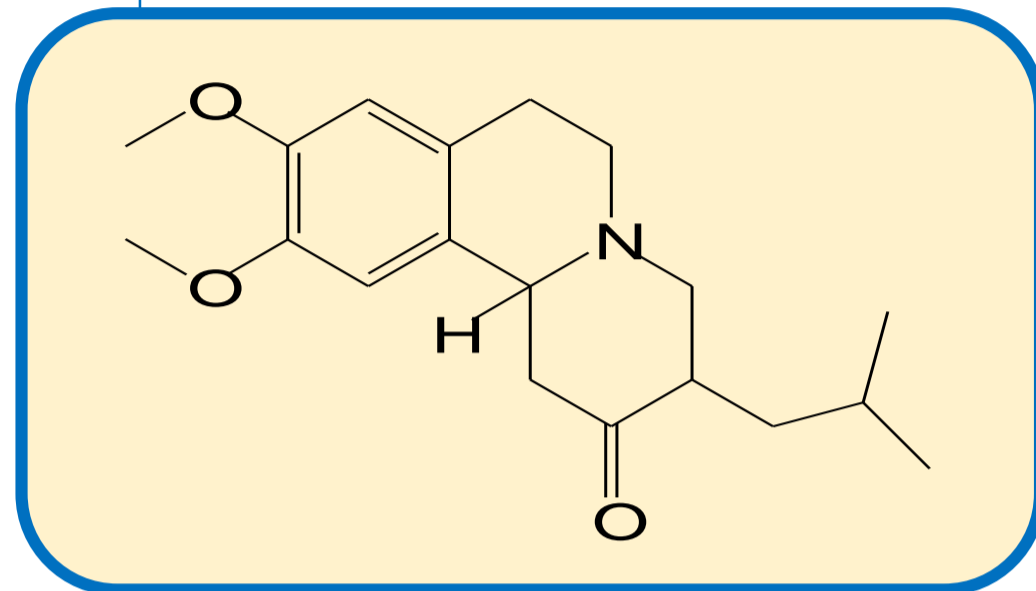

SEPARATION PERFORMANCE

t_{R1} : 5.2	N_1 : 9560	T_1 : 1.2	
t_{R2} : 7.1	N_2 : 8456	T_2 : 1.3	R_s : 7.2

