

Sotalol

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87325

Mobile phase: n-Hexane : Ethanol : DEA (50/50/0.1, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 230 nm

Diluent: Ethanol

Sample Concentration: 1.0 mg/mL

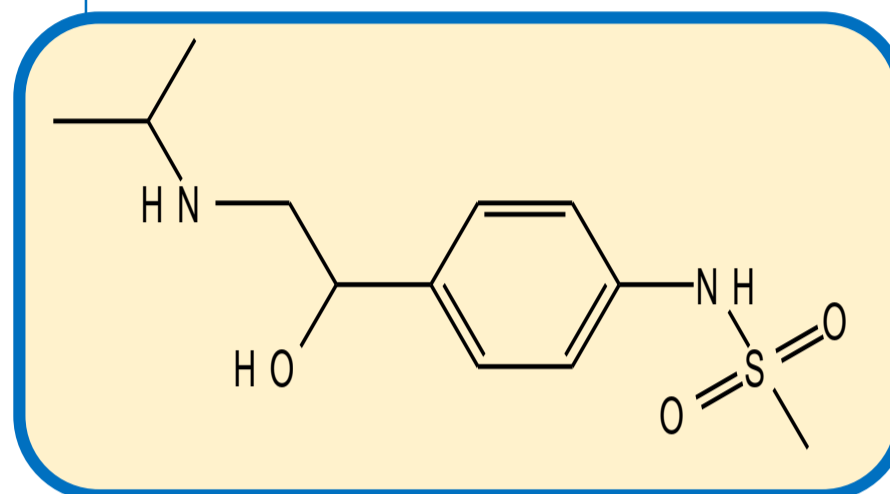

SEPARATION PERFORMANCE

t_{R1} : 4.2	N_1 : 7152	T_1 : 1.6	
t_{R2} : 5.3	N_2 : 8323	T_2 : 1.5	R_s : 4.9

