


DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Clobenzorex

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 150mm, 5 µm]
Part No: 87324
Mobile Phase: 10mM Ammonium bicarbonate (aq.) : ACN (50/50, v/v)
Column Oven Temperature: 25°C
UV Detection: 210 nm
Flow Rate: 1.0 ml/min
Diluent: Water/ACN (1:1)
Sample Concentration: 2.0 mg/ml

Clobenzorex

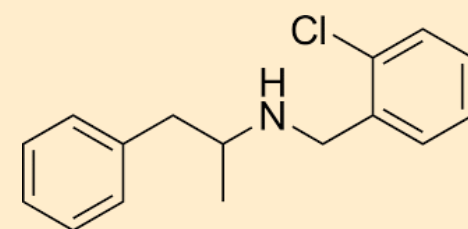

SEPARATION PERFORMANCE

t_{R1} : 11.7	N_1 : 4536	T_1 : 1.0	
t_{R2} : 13.6	N_2 : 4478	T_2 : 1.1	R_s : 2.6

